

Operation Manual

RyfeMax Platinum



Table Of Contents

Table of Contents

Introduction:	4
Basics:	5
What is electrotherapy treatment for?	5
How does a Rife machine work?	5
Technical Specifications:	6
How to use the RyfeMax Platinum machine:	8
How to connect the power?	8
Quick Start:	9
Operational flow using the RyfeMax Platinum machine:	9
Add a new client:	10
1. Click "Client details"	10
2. Click "Add a new client"	10
3. Type in the client details:	10
Click on each field you wish to add details for:	10
Type in the details:	10
Click "Accept"	10
4. Click "Save"	10
Click a field:	10
Add an ailment to favourites:	11
Running all of your favourites:	12
Running a single favourite:	13
How to use the RyfeMax Platinum 100% contact-less device:	14
Note – To operate the touch screen:	14
Resume feature:	14
Power save mode:	14
How To Get The Best Out Of Your RyfeMax Platinum machine:	15
Flowchart:	15
How often do I need to use the RyfeMax machine?	16
Settings menu:	17
How do I adjust the sleep duration:	17
How do I adjust the duration?	18
How do I lock the screen?	19
How to unlock the screen:	19
How do make the ailments I am running repeat?	20
Run multiple ailments menu:	21
Sweep All HC Parasites (Dr Hulda Clarke parasite frequencies):	21
Sweep all Cancers:	22
Sweep all Cold & Flu:	23
Sweep all Corona Virus COVID-19:	24
Sweep all favourites:	25
Sweep all known frequencies:	26
Sweep all parasites:	27
Sending ailments to the RyfeMax Mini and Mini+:	28
Connecting the cables:	28
Viewing ailments on the RyfeMax Mini and Mini+:	29
Special considerations:	30
Selecting the correct Rife frequencies:	30
What is the herxheimer reaction?	30
What to do when you get a herxheimer reaction?	30
Warnings:	30

Suggestions:.....	31
General:.....	31
For best results:.....	31
Allergies:.....	32
Arthritis Rheumatoid:.....	32
Arteriosclerosis:.....	32
Back pain:.....	32
Cancer:.....	33
Diet for People with cancer and diabetes:.....	34
Eating meat:.....	34
Bicarbonate of soda:.....	34
Candidiasis (candida overgrowth, yeast infection, thrush):.....	35
Cold Sores:.....	35
COVID-19:.....	35
Dementia:.....	35
Detox:.....	36
Diabetes type 1 (or late onset type 1):.....	36
Diabetes type 2 (or late onset type 1):.....	37
Epstein-barr Virus infections:.....	37
Fatty liver:.....	37
Heart diseases:.....	37
High blood pressure:.....	38
Osteoarthritis:.....	38
Prostate enlarged:.....	38
Myocarditis:.....	38
Parasites:.....	38
Shingles:.....	39
Sleep problems – getting to sleep:.....	39
Thrombosis.....	39
Thrombocytopenia.....	39
Thyroid low, Hashimoto's or reverse T3.....	40
Tinnitus:.....	40
Troubleshooting.....	41
DIAGNOSTIC FAILURE WINDOW is displayed.....	41
I have a program running but the buttons do not work.....	41
I can't find my ailment on the list.....	41
The machine seems to make a noise sometimes.....	41
Sales And Enquiries:.....	41

Introduction:

Congratulations on the purchase of your RyfeMax Platinum machine. The following guide provides an introduction to Electro-therapy and recommendations for the use of your device. We wish you every success with your therapeutic treatments and the positive experience of using your new device!

Only use the RyfeMax Platinum machine in accordance with the instruction guide.

Ensure that water and other fluids do not come in contact with the RyfeMax Platinum machine.

DISCLAIMER:

RyfeMax and any of its associated products do not diagnose, cure, prevent or treat disease. It is an experimental device.

If you have a medical condition or concern please consult an appropriate health care professional.

By using this device you are accepting these conditions and in no way hold RyfeMax and its affiliates responsible for any outcomes.

All suggestions in this document relating to using the RyfeMax machine for diseases or medical conditions, should in no way be interpreted as diagnosing, treating or preventing such disease or medical condition.

Basics:

What is electrotherapy treatment for?

Bacteria, viruses and parasites cause most of today's diseases.

Electro-therapy is used to kill bacteria, viruses and parasites. There are also various frequencies that have been found to be beneficial in aiding the body's repair mechanisms.

Electrotherapy is also used for many other conditions that are not caused by bacteria, viruses or parasites.

How does a Rife machine work?

Dr Royal Rife discovered that all organisms (in fact all things) have a resonant frequency.

An example of this was the collapse of the Broughton suspension bridge that collapsed in 1831 caused by marching soldiers whose footsteps matched exactly the resonant frequency of the bridge.

Another example, is in the medical treatment of Kidney stones. Generally someone with Kidney stones goes to hospital and is placed in a bath. The patient is then subjected to ultrasonic frequencies that match the resonant frequency of the kidney stones and they explode into gravel, which is then passed during urination.

Dr Rife discovered, when organisms such as parasites, bacteria and viruses are subjected to their unique resonant frequency, they explode! Dr Rife (and others) discovered the organisms that are the underlying cause of a disease and the resonant frequency that caused the organism to be destroyed.

Further research has been carried out by a number of other scientists who have discovered the resonant frequency of many more organisms that are the underlying cause of many of our health issues. Resonant frequencies for over 3,500 diseases have been discovered so far.

As we hold a copy of Dr Rife's original research notes, clinical trials, double blind studies, peer reviews, patient testimonials and circuit diagrams, our RyfeMax Platinum machine is based on Dr Rife's work and as such is an accurate representation of his device, using modern day electronics.

Technical Specifications:

FREQUENCY RANGE

1 hz - 4,000,000 hz (4 megahertz)

ACCURACY

1 part in 20 million

DESIGN

The RyfeMax Platinum machine uses a state of the art digital microprocessor, 10 x analog processors and a touch screen.

PRE-PROGRAMMED AILMENTS

The Rife frequencies for over 3,500 ailments have been pre-programmed into the RyfeMax Platinum machine. No need to add the frequencies in yourself, just search for the ailment you are interested in, and click Run.

FAVOURITES

The number of favourites is only limited to the amount of memory storage available, and is in excess of 100,000.

VOLTAGE INPUT

5v DC 5 amps from a USB PD-3.0 100w or more power device. Make sure that you ONLY use a USB PD-3.0 100w or more charger. A suitable USB charger is included in the retail pack.

There are many non-genuine, 3rd party USB phone chargers available on the market, some even claim 1 amp or more output, but in fact output less than 1/10th of the 1 amp power requirement.

If the Screen of the RyfeMax Platinum continually re-boots or glitches, then the USB charger is most likely **NOT** suitable.

OUTPUT FREQUENCY

The RyfeMax Platinum device outputs a radio frequency carrier wave that is modulated with an ailment specific square wave.

RIFE FREQUENCIES

The frequencies used by the RyfeMax Platinum machine are from the ETDFL Association.

ETDFL is a project based in Berlin Germany, with a vision to contribute innovative information and technology to the world.

ETDFL Have been compiling disease frequency data lists from 12 Bio-resonance Clinics across the planet since June 2006.

Data is sourced from 12 Clinics which use the Quantum SCIO Bio-resonance machine to record Disease frequencies of actual clients, and then produce the 10 most accurate frequencies which are deemed suitable for treatment with the Rife Machine, Frequency Generator Devices, or

The ETDFL project began in 2006 when the original CAFL frequency list was discontinued. ETDFL became the industry standard for Bio-resonance Frequencies and is now used by more than 100,000 Rife Machine enthusiasts and clinics across the planet. The accuracy of the frequencies produces positive results for the clients who use this Frequency List in clinics and homes in every country.

ETDFL is updated once a year. Latest PDF versions are available at this website. The 12 International Language PDF versions of ETDFL will allow you to search the ETDFL Database for Disease frequencies which may be more up to date than the Paperback ETDFL provided with your Rife Machine.

SAFETY

There are 3 levels of safety built into the machine, including electronic circuitry and monitoring of the output. This ensures that the output voltage remains within manufacturers specifications.

DO NOT USE THE MACHINE IN OR ON YOUR BED, ON OR UNDER ANY PILLOWS, BLANKETS OR OTHER ITEMS THAT MAY CAUSE THE UNIT TO OVERHEAT.

How to use the RyfeMax Platinum machine:

How to connect the power?



- ✓ Make sure the power display of the USB-C cable is connected to the RyfeMax Platinum machine and is facing upwards.
If the USB-C cable is connected upside-down, power will not be delivered to the RyfeMax Platinum machine.

Quick Start

Operational flow using the RyfeMax Platinum machine

1. Add a new client

- ✓ A default client named A Business has already been set up for you
- ✓ If the client you wish to use has already been set up, you can omit this step



Client details

2. Add ailments for the selected client (or optionally run a single ailment)



Select and run favourite ailments

3. Run all of the favourites for your selected client



Run multiple ailments

Add a new client

1. Click "Client details"
2. Click "Add a new client"
3. Type in the client details
 - ✓ Click on each field you wish to add details for
 - ✓ Type in the details
 - ✓ Click "Accept"
4. Click "Save"



Click a field



Add an ailment to favourites

2. Click "Select and run favourite ailments"
3. If you need to select the client, then Select the client
 - ✓ Click "Select client"
 - ✓ Select the client



Select client

4. Click "Add new ailment"
5. Type in the first few letters of the ailment you wish to add
 - ✓ For a list of ailments, go to www.ryfemax.au and click the "Ailments" button
6. Click Accept
7. Select the ailment you wish to add
8. Click Save
9. Repeat steps 4 – 8 until all the ailments you need for this client are added

5. Click "Done"



Select and run favourite ailments



Add new ailment

Type first few letters of the ailment
Select the ailment

Repeat the above steps until all ailments have been added



Done

Running all of your favourites

1. If you need to select the client, then Select the client

- ✓ Click "Select and run favourite ailments"
- ✓ Click "Select client"
- ✓ Select the client
- ✓ Click "Done"



Select and run favourite ailments



Select client



Done

2. Click "Run multiple ailments"

3. Click "Sweep all favourites"

4. Click "Run"



Run multiple ailments

Sweep all favourites



Run

Running a single favourite

1. Click "Select and run your favourite ailments"
2. If you need to select the client, then Select the client
 - ✓ Click "Select and run favourite ailments"
 - ✓ Click "Select client"
 - ✓ Select the client
 - ✓ Click "Done"



Select client

3. Select the favourite you wish to run

4. Click "Run"



Select and run favourite ailments

Select the favourite you wish to run



Run

How to use the RyfeMax Platinum 100% contact-less device

The RyfeMax Air 100% contact-less machine works within a 10 metre radius, line of sight. This being said, the closer you are to the device, the stronger the signal.

Each 1 metre you move away from your RyfeMax, the signal strength is degraded by 50%.

Note – To operate the touch screen:

To operate the touch screen, use either your finger or a rubber pointing device.

Resume feature

When moving the machine from location to location, there is no need to click any buttons or do anything on the screen.

1 minute after powering on the RyfeMax machine, it will automatically resume at the beginning of your most recent 10 min frequency block.

Power save mode


When running a program or ailment, after a few seconds the screen will go blank. This is normal operation as it increases the longevity of the screen.

To refresh the screen, simply touch the screen.

How To Get The Best Out Of Your RyfeMax Platinum machine


Flowchart:

How to get the best out of your RyfeMax machine

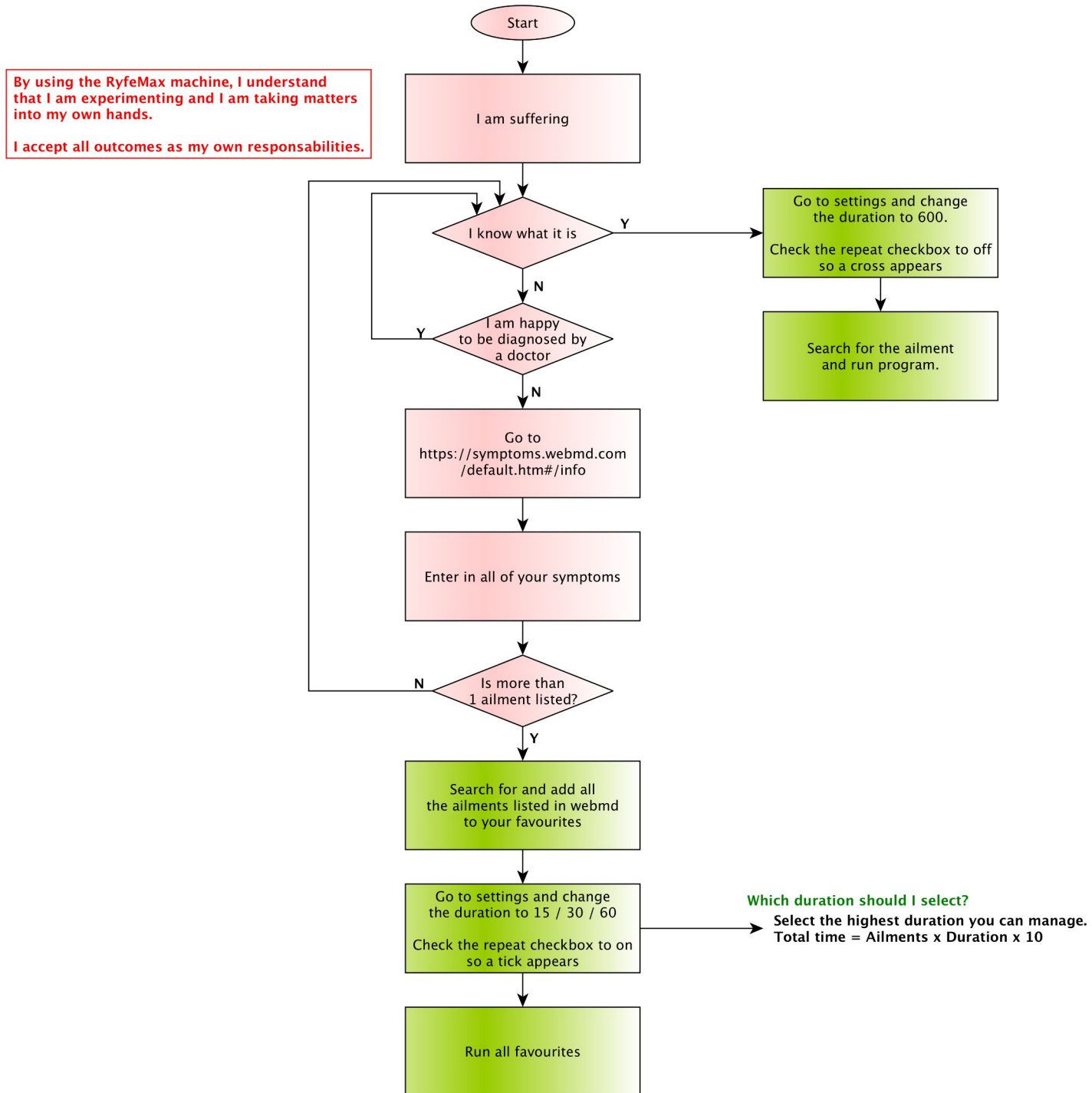

How often do I need to use the RyfeMax machine?

How Often?

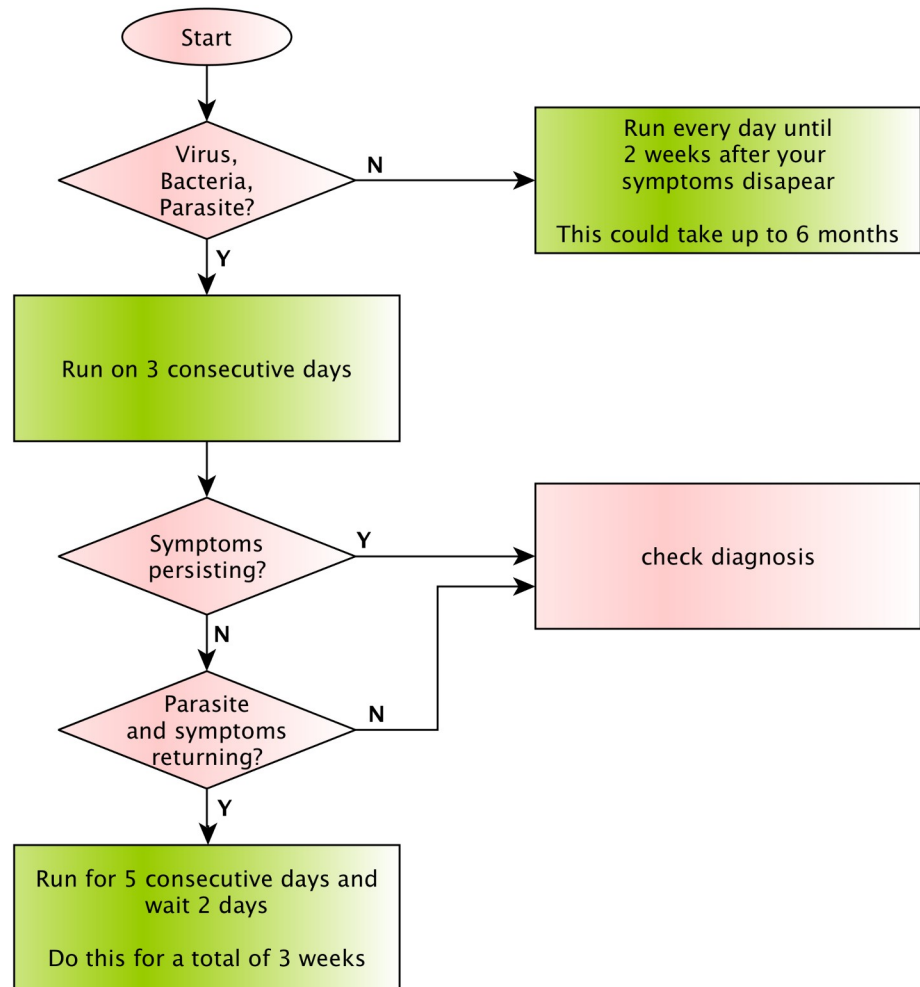

Settings menu

How do I adjust the sleep duration

1. Click "Settings"
2. Click the "Sleep" field (The box next to Sleep)
3. Type in the amount of time before the screen sleeps in seconds
4. Click "Accept"
5. Click "Save"



Click the Sleep field



How do I adjust the duration?

Usually it is best to set the duration field to 600

1. Click “Settings”
2. Click the “Duration” field (The box next to Duration)
3. Type in the duration of time each frequency should run in seconds
4. Click “Accept”
5. Click “Save”



Click the Duration field



How do I lock the screen?

1. Click "Settings"
2. Click the "Auto Lock" field (The box next to Auto Lock) so it becomes green
3. Click "Save"



Click the Auto Lock field



How to unlock the screen

1. Power off
2. Wait 10 seconds
3. Power on
4. Click the logo on the screen
5. Click "Settings"
6. Click the "Auto Lock" field (the box next to Auto lock) so it becomes clear
7. Click "Save"



Click the Auto Lock field



How do make the ailments I am running repeat?

1. Click "Settings"
2. Click the "Repeat" field (The box next to Auto Lock) so it becomes green
3. Click "Save"



Click the Repeat field



Run multiple ailments menu

This folder contains a number of options that are composites of any one ailment's frequencies.

Sweep All HC Parasites (Dr Hulda Clarke parasite frequencies)

Dr Hulda Clarke in her book "The cure for all cancers" listed a number of parasites that transport the cancer viruses into your body.

This program delivers all the frequencies for these parasites.

1. Click "Run multiple ailments"
2. Click "Sweep All HC Parasites"
3. Click "Run"



Sweep All HC Parasites



Sweep all Cancers

It is advisable that you follow the instructions set out in the Special Instructions Cancer section.

This program is only if you know you have Cancer, but do not know which Cancer you have.

This selection will go through all the known Rife frequencies for the different cancers. There are 36 different cancer frequencies, so take this into consideration, when running sweep all cancers.

Continue to do this until the entire sweep of Cancer frequencies is complete.

1. Click "Run multiple ailments"

2. Click "Sweep All Cancers"

3. Click "Run"



Sweep All Cancers



Sweep all Cold & Flu

This will go through Common Cold, Influenza A, Influenza B, Influenza, Bird Flu, Swine Flu.

You may find you are bringing up green phlegm for the next few days.

1. Click "Run multiple ailments"
2. Click "Sweep All Cold And Flu"
3. Click "Run"



Sweep All Cold & Flu



Sweep all Corona Virus COVID-19

This will go through the frequencies for COVID.

1. Click "Run multiple ailments"
2. Click "Sweep All Corona virus COVID-19"
3. Click "Run"



Run multiple ailments

Sweep All Corona virus COVID-19



Sweep all favourites


This will go through all the ailments in your favourites list, and pulse each frequency of your favourite for a period set in the “Duration” value of settings.

When you run “Sweep All Favourites”, the RyfeMax Platinum machine will run every frequency for the first ailment on your favourites list one after another, then do the same for the next ailment on your favourites list and so on until all of the frequencies for all of the ailments on your ailments list have been run.

If you have the Repeat checkbox ticked in Settings, then this process will continue on indefinitely.

Running “Sweep All Favourites” can take some time depending on the number of ailments you have set up in your favourites (minus the “Sweep All programs” and “Earth Resonant Frequency”) and the duration you have set up in Settings.


1. Click “Run multiple ailments”
2. Click “Sweep All Favourites”
3. Click “Run”

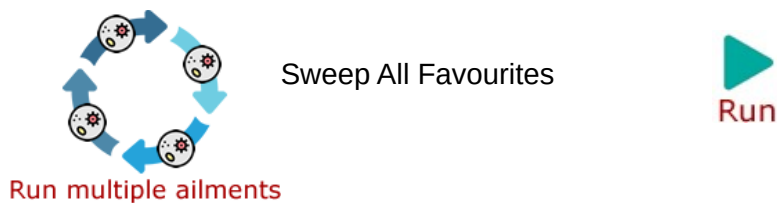

Sweep all known frequencies

This selection will go through all the known beneficial Rife frequencies.

1. Click “Run multiple ailments”
2. Click “Sweep All Known Frequencies”
3. Click “Run”

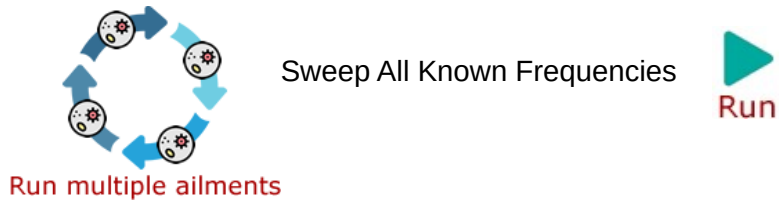

Sweep all parasites

This selection will go through all the known Rife frequencies for the different parasites. There are 390 different parasite frequencies, so take this into consideration, when running sweep all parasites.

When you are finished for the time being, simply power off the RyfeMax Platinum machine while it is running. The next time you power up the RyfeMax Platinum machine, it will resume at the start of the most recent 10 minute block where you left off.

Continue to do this until the entire sweep of Parasite frequencies is complete.

1. Click "Run multiple ailments"

2. Click "Sweep All Parasites"

3. Click "Run"



4. Sweep All Parasites



Sending ailments to the RyfeMax Mini and Mini+

Your RyfeMax Platinum machine is equipped with the ability to also send your favourite ailment frequencies to the RyfeMax mini.

1. Connect power to both the RyfeMax Platinum and the RyfeMax Mini
2. Connect the gold end of the supplied cable to the Ryfemax Platinum
3. Connect the other end of the supplied cable to the RyfeMax Mini
4. Click “Export ailments”
5. Select what you wish to export



Select what you wish to export

Connecting the cables



Viewing ailments on the RyfeMax Mini and Mini+

Your RyfeMax Platinum machine is equipped with the ability to also view which ailments have been set up on the RyfeMax Mini and RyfeMax mini+.

This will only work if the ailments were loaded on the RyfeMax Mini or RyfeMax Mini+ using the RyfeMax Platinum machine.

1. Connect power to both the RyfeMax Platinum and the RyfeMax Mini
2. Connect the gold end of the supplied cable to the Ryfemax Platinum
3. Connect the other end of the supplied cable to the RyfeMax Mini
4. Click "Export ailments"
5. Select "View Mini favourites"



Select View Mini favourites

Special considerations

Selecting the correct Rife frequencies

If you are unsure of your condition, then we recommend visiting your doctor or medical professional to get a diagnosis.

If you use the wrong Rife frequencies by selecting the wrong ailment, nothing will happen. Of course the ailment you wish to use the RyfeMax Platinum machine machine for will not be addressed

You cannot “do anything wrong” with using the RyfeMax Platinum machine. If (for example) you use the device to sort out a particular parasite, and if you do not have this parasite, then the electrical frequency may simply fade without effect on the body.

What is the herxheimer reaction?

When you begin using the RyfeMax Platinum machine, you may experience the “Herxheimer effect” or “Die off syndrome”, which manifest themselves as nausea, fatigue, dizziness and flu-like symptoms.

This brief overloading of the liver and other eliminating organs worsen your symptoms at first. This can be a for a few hours, but can also last up to a few days before you begin to feel much better.

It is important to know that the worsening of symptoms does not mean what you are running on the RyfeMax Platinum machine has failed; in fact, it is exactly the opposite, namely the evidence for a very strong and effective detoxification.

What to do when you get a herxheimer reaction?

The most important thing is to accept this reaction, as a confirming result of what you are running on the RyfeMax Platinum machine. If it is very severe, then reduce the time using the RyfeMax Platinum machine, drink lots of water (at least 3 L per day) and rest (no exercise!), eat lightly to not burden your digestive system additionally, add freshly squeezed lemon juice to your drinking water and make sure that your digestion system is working properly.

Warnings

If you use the frequencies of Lactobacillus Acidophilus then you may have to also take a pro-biotic.

Do not place the machine on a soft surface such as a pillow or carpet as doing so may cause the machine to overheat and present a fire risk.

Suggestions:

In this chapter, we will list some suggestions on how to get the best results using your RyfeMax Platinum machine. These suggestions are by no means exhaustive, they are merely some of the more common ailments people suffer from.

Our website is www.ryfemax.au

Pay particular attention to the flowchart on pages 15 and 16. This helps people to understand how to manage ailments.

For a list of ailments, go to <https://www.ryfemax.au/DisplayPage?&Name=Ailment>

If you can't find the ailment on the list, then lookup google for other names for the ailment, then search again on our website.

General:

If you are wondering if the machine is working or not, the machine has 3 levels of output checking and will alert you with an error window if there is a problem.

If you get no improvement in your symptoms, then it means the diagnosis of your problem is incorrect.

If you are suffering from a migraine and you use the machine for brain cancer nothing will happen, good or bad. I.E. you need a correct diagnosis!

Case in point. Leone and I contracted Malaria in Ecuador. We had really bad flu symptoms and used the machine for cold & flu, gram positive, gram negative, influenza A, Influenza B, nothing seemed to work.

We went to www.webmd.com and entered our symptoms. It came up with only 1 thing Malaria. 1st time running Malaria the filthy headache was gone, after the 3rd time running Malaria all symptoms were gone minus the fatigue, which went away after our bodies regenerated enough red blood cells which took about 3 weeks.

You can not use the machine too much! The more time on the machine the quicker the results.

For best results:

Run the machine daily with the duration set to 600. (doesn't have to be at the same time of day).

If you have an ailment like chronic fatigue, then run the machine on that ailment daily.

Also add possible root causes for the ailment, such as Ross River fever (alpha virus infections), dengue, Epstein Barr or Malaria. Look into your history as to any of those kinds of viruses that "stay with you" through your life and add them to favourites.

Allergies:

1. Hypersensitivity Immediate
2. Hypersensitivity
3. Rhinitus

Arthritis Rheumatoid:

1. Epstein-Barr Virus Infections
2. Arthritis Rheumatoid

Arteriosclerosis:

1. Arteriosclerosis

If you are on blood pressure medication, then you MUST monitor your blood pressure every day, and if necessary adjust your medication.

Back pain:

Back pain can have a number of causes, it's a matter of trying them and see which one improves your symptoms. Run either of the ones listed below.

1. Spondylolithesis
2. Intervertebral Disk Displacement
3. Nerve Compression Syndrome
4. Back Pain

You should see improvement within a day or two. If not, then try one of the other options.

Cancer:

Cancer is a war being played out in your body. Scientists have discovered that Cancer feeds on Sugar and sugar substitutes and grows prolifically in an acid environment, but dies off in an alkaline environment. Keeping these things in mind, we suggest the following...

1. Bx and By Virus
2. Cancer ... (where ... is the name of the primary cancer)
3. Carcinoma
4. If it is a Squamous cell carcinoma then Cancer Metastasized Squamous Neck
5. Run these programs 24/7 until 1 month after all the symptoms are gone

In some cases it might be better to use a dedicated portable Rife machine such as the Ryfemax Air professional, RyfeMax Air standard, RyfeMa Mini or RyfeMax Mini+ models so you can run the programs 24/7

Also...

1. Follow a strict diet which removes all sugars, fake sugars, take out and seed oils from your diet.
2. Add at least 200gm of beef, lamb, venison, goat or rabbit every day before 3pm
3. Eat all your carbs (without protein) after 3pm
4. Eat nothing after 6pm
5. Take bi-carbonate of soda (see instructions below)

Diet for People with cancer and diabetes:

No:====> pasta, rice, bread, flour, junk food, sugar or sugar substitutes or anything containing dextrose, sugar, high fructose corn syrup, soft drinks (including the zero sugar varieties), coffee, alcohol, biscuits, chips, avoid all processed foods, so nothing from a box, can or packet.

Avoid all seed and vegetable oils. Use real butter (the one in the paper, not blends in tubs), ghee or coconut oil.

No take out food.

Basically eat a diet of meat and veg (spices are ok) or meat and salad. The preference is organic if this is possible but not mandatory.

Whole fruit is OK. Keep away from Genetically Modified foods such as corn.

Honey is OK (must be raw honey from a farmer, not the supermarkets fake honey which is processed High Fructose Corn Syrup and is labelled Honey, it comes from China) Other supermarket Honey is heat treated which changes it's chemical composition.

Eggs, milk, cheese and non wheat flour is ok.

If you have sugar cravings, then take Chromium, but the best time is to take it at 2pm daily.

Eating meat:

Eat as much meat (beef, lamb, venison, goat or rabbit) as you can (not processed meats like ham or salami. and DO NOT eat pork or chicken).

When eating meat, either add the juice of 2 freshly squeezed lemons (not the fake lemon juice from the supermarket) to your meal or sprinkle 1/2 teaspoon Celtic Sea Salt on your meat. These actions will ensure that your body maintains an alkaline state.

ONLY USE CELTIC SEA SALT - no substitutes.

ABSOLUTELY NO CHEATING:

Bicarbonate of soda:

Take 1 teaspoon of Bicarbonate of soda in water 2 times a day. 3 weeks on, 2 weeks off.

Taking the Bicarbonate of soda 3 weeks on and 2 weeks off, greatly reduces the impact on your digestive system, and is critical in dealing with Cancer.

Monitor your p.h. every day using the PH strips that go under your tongue.

If your p.h. goes above 8.2 then reduce the amount of Bicarbonate of soda until it goes below 8.0.

Candidiasis (candida overgrowth, yeast infection, thrush):

1. Candidiasis

Usually need to run Candidiasis for 30 days consecutive, duration = 600. We have found this generally causes really bad flu like symptoms because of the toxins released via the Candida die off.

If the symptoms become unbearable, then run Candidiasis one day, Liver Detox the next and have a rest the third day. Of course, it will then take longer.

While running Candidiasis, we suggest taking MRM Cardio Chelate from iHerb. Take 2 tablets 3 times a day (or whatever your healthcare professional suggests). This will bind to the toxins and facilitate elimination.

If you look at your tongue, if it is no longer white, then your Candida overgrowth has subsided. Run the Candidiasis for another 15 times.

Cold Sores:

1. Herpes labialis

Running programs for certain viruses such as yellow fever, shingles, cold sores, dengue will activate the virus then kill it, so you will experience some symptoms while using the machine.

COVID-19:

If you have the disease, the run "Sweep all corona virus COVID-19" from the "Run multiple ailments menu".

As a preventative, run Coronaviridae Infections

If you have been subjected to others shedding the ACE-2 Spike protein, run Coronaviridae Infections

DEMENTIA:

1. Dementia

2. Dementia Vascular

Detox:

When running any sort of Detox, you **MUST** also take a chelating agent such as CH-77 or MRM Cardio Chelate from iHerb. Take 2 tablets 3 times a day (or whatever your health care professional suggests) while detoxing.

If you start to get a rash while detoxing, this means that your body is eliminating the toxins via your skin. If necessary, do a detox one day, then have 2 days break (but still take the MRM Cardio Chelate).

Diabetes type 1 (or late onset type 1):

All diabetics have a common fluke parasite, *Eurytrema pancreaticum*, the pancreatic fluke of cattle, in their own pancreas. It seems likely that we get it from cattle, repeatedly, by eating their meat or dairy products in a raw state. It is not hard to kill with a machine but because of its ineffective stages in our food supply we can immediately be reinfected.

Eurytrema will not settle and multiply in our pancreas without the presence of wood alcohol (methanol). Methanol pollution pervades our food supply—it is found in processed food including bottled water, artificial sweetener, soda pop, baby formula and powdered drinks of all kinds including health food varieties.

I presume wood alcohol is used to wash equipment used in manufacturing.

If your child has diabetes, use nothing out of a can, package or bottle except regular milk, and no processed foods.

1. Diabetes Melitis
2. Parasites *Eurytrema pancreaticum*
3. Create a custom ailment with 421.000 in each cell
4. Strictly maintain the diet under the section “Diet for people with cancer and diabetes”

Diabetes type 2 (or late onset type 1):

All diabetics have a common fluke parasite, Eurytrema pancreaticum, the pancreatic fluke of cattle, in their own pancreas. It seems likely that we get it from cattle, repeatedly, by eating their meat or dairy products in a raw state. It is not hard to kill with a machine but because of its ineffective stages in our food supply we can immediately be reinfected.

Eurytrema will not settle and multiply in our pancreas without the presence of wood alcohol (methanol). Methanol pollution pervades our food supply—it is found in processed food including bottled water, artificial sweetener, soda pop, baby formula and powdered drinks of all kinds including health food varieties.

I presume wood alcohol is used to wash equipment used in manufacturing. If your child has diabetes, use nothing out of a can, package or bottle except regular milk, and no processed foods.

1. Diabetes Melitis
2. Parasites Eurytrema pancreaticum
3. Create a custom ailment with 465.000 in each cell
4. Strictly maintain the diet under the section "Diet for people with cancer and diabetes"

Epstein-barr Virus infections:

The Epstein-barr virus has over 60 variations which are the root cause of many, many ailments.

1. Epstein-barr virus infections

Fatty liver:

1. Fatty liver

Heart diseases:

1. Heart diseases

If you are on blood pressure medication, then you MUST monitor your blood pressure every day, and adjust your medication if your blood pressure goes below 120 / 80.

High blood pressure:

1. Hypertension
2. Hypertension portal
3. Hypertension pulmonary
4. Heart diseases
5. Low potassium

If you are on blood pressure medication, then you **MUST** monitor your blood pressure every day, and your medication if your blood pressure goes below 120 / 80.

Osteoarthritis:

1. Arthritis
2. Critical Boron and critical Magnesium deficiency (see Dr Walter Last)

The Arthritis program covers Osteoarthritis.

Prostate enlarged:

1. Prostate enlarged
2. Prostate infection / pain

Myocarditis:

1. Myocarditis
2. If you have been vaxxed for Covid-19, also use Coronaviridae Infections

Parasites:

According to Wikipedia, 30% - 50% of people have parasites (yes even in the western world), other sites such as <https://www.victoriahealth.com/editorial/most-of-us-have-intestinal-parasites> quote as high as 80% of people have parasites.

We like to run Parasites General Comprehensive for 3 consecutive days before doing anything else.

Generally run the specific parasite, or Parasites General Comprehensive for 3 consecutive days duration = 600.

Shingles:

1. Herpes Zoster

If the symptoms do not go away, then either you are immune compromised, you have other herpes viruses in the herpes family or the Herpes Zoster virus has damaged the nerves.

If you are immune compromised or you have other herpes viruses in the herpes family, run Herpesviridae Infections and Herpes Zoster every day until the symptoms disappear (this could take many weeks).

If you have nerve damage, then run Neuralgia. This could take many months, as the repair of nerve is a very slow process.

Running programs for certain viruses such as yellow fever, shingles, cold sores, dengue will activate the virus then kill it, so you will experience some symptoms while using the machine.

Sleep problems – getting to sleep:

1. Earth Resonant Frequency.

Run this when you are trying to get to sleep.

Thrombosis

1. Thrombosis

2. 1/8 Teaspoon Cayenne pepper directly on the tongue 1t thing in the morning

If you have been vaxxed for COVID-19

1. Coronaviridae Infections

Thrombocytopenia

1. Thrombocytopenia

If you have been vaxxed for COVID-19

1. Coronaviridae Infections

Thyroid low, Hashimoto's or reverse T3

1. Epstein-barr virus infections

2. Hypothyroidism

3. Adrenal gland diseases

4. Combat disorders

This will take a year or more

Tinnitus:

1. Tinnitus

2. Hearing disorders

3. Epstein-barr virus infections

Troubleshooting

DIAGNOSTIC FAILURE WINDOW is displayed

The RyfeMax Platinum machine has 3 levels of output checking to ensure that the machine is working correctly.

However, if a problem has arisen in the output, then this message will show up in red.

1. Try powering the machine off, wait 10 seconds and power the machine back on again.
2. If the above does not work, then your machine will have to be returned to the factory for repair.

The turn around time is usually 1 -2 working days plus time for postage both ways.

I have a program running but the buttons do not work

If in settings, you check the Autolock feature on, then this will disable the Main button when running a program. To gain access to your machine again, do the following..

1. Power the machine off
2. Wait 10 seconds
3. Power back on.
4. At the disclaimer window, click the Logo
5. You now have access to the favourites window

I can't find my ailment on the list

Sometimes, the ailment you are interested in is spelled a little differently or goes by another name.

If you can't find your ailment on the the list, then search on Google for other names of the ailment and search on our website <https://www.ryfemax.au/DisplayPage?Name=Ailment> for the other names.

The machine seems to make a noise sometimes

Sometimes, the machine will make a noise. This is only for certain frequencies and when the screen is lit up and is perfectly normal.

Sales And Enquiries:

We are always looking for distributors. For sales and enquiries please contact:

info@ryfemax.com